

Lab Scale Process Vessels

Three different vessel sizes, with low hold-up volume - the key to minimizing overall system holdup. All have a conical bottom and a low point drain.

Small Scale Vessels


- Designed with a conical base and low-point drain that minimizes liquid hold-up
- Female luer lock outlet at base; this facilitates easy integration to any process flow path
- Lid with 3 holes easily inserts into vessel - 2 holes for placement of dip tubes and one for venting
- Locator at flat bottom for stir bar for process mixing
- Transparent with graduation marks for accurate process volume measurement
- Made of polysulfone which is compatible with many chemicals including sodium hydroxide
- Available in 2 different sizes of 600mL and 140mL

Contact Materials
 Vessel: Polysulfone
 Dip Tubes: Nylon
 Cap: TPV



600mL with 40mL in conical section

140mL with 12mL in conical section

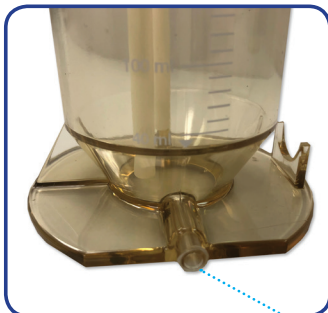


Luer Outlet




Stir bar & polysulfone 3-way stopcock included

Dip tubes that can easily be positioned as needed

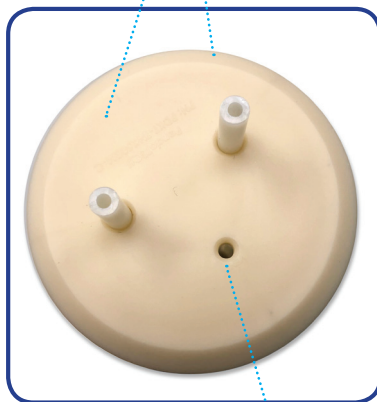

Vent hole sized as a female luer fitting in case a vent filter is desired to be used

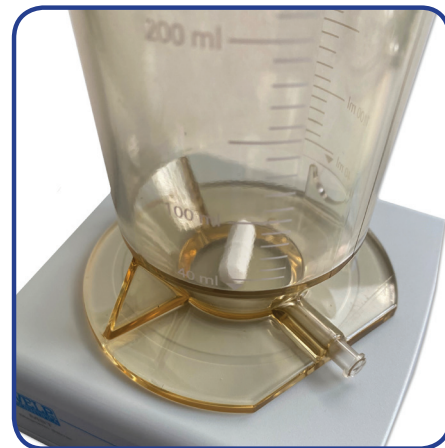


Pocket for the included magnetic stir bar*

*Stir plate sold separately

PendoTECH® is a registered trademark of PendoTECH, all rights reserved.

2 Liter Vessel

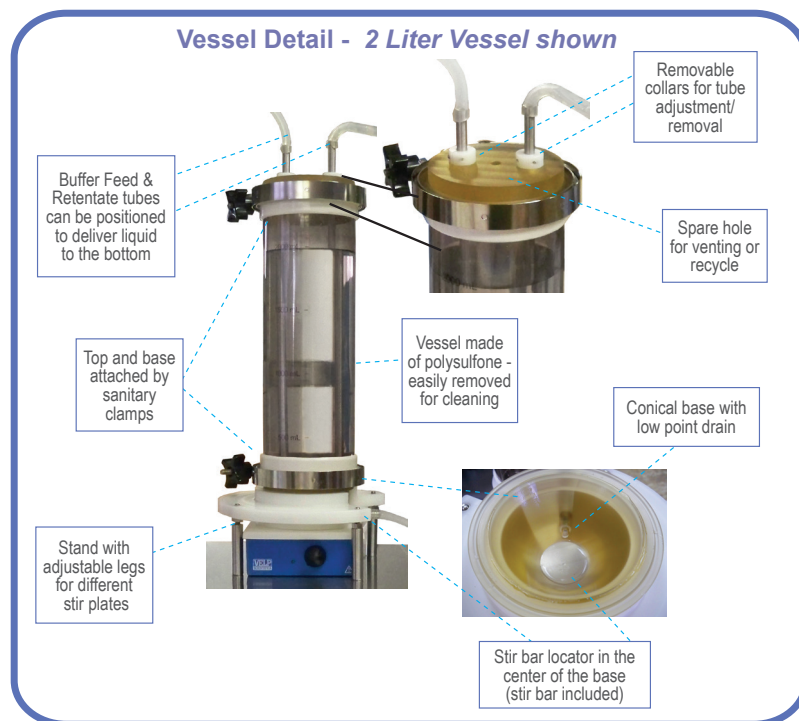

Ordering Information

PDKT-TNK	Complete 2L tank setup with stir plate, 2 dip tubes, and a stand
PDKT-TNK500M	600mL - conical bottom vessel (polysulfone) with low point drain, includes lid, 2 dip tubes, 1 polysulfone 3 - way stopcock, and stir bar - (DOES NOT INCLUDE STIR PLATE)
PDKT-TNK125M	140 mL - conical bottom vessel (polysulfone) with low point drain, includes lid, 2 dip tubes, 1 polysulfone 3 - way stopcock, and stir bar - (DOES NOT INCLUDE STIR PLATE)
TNK-2L-STRPLATE	Stirrer plate

Warranty/Disclaimer

All descriptions, representations and/or other information concerning Product on the PendoTECH website and/or contained in PendoTECH's advertisements, brochures, promotional material, or statements made by employees or sales representatives of PendoTECH are solely for general informational purposes only and are not binding upon PendoTECH. No employee or sales representative of PendoTECH shall have any authority to establish, expand or otherwise modify PendoTECH's warranty associated with the sale of Product. PendoTECH, LLC makes no warranty or representation regarding whether or not a customer's end use of any PendoTECH product or system infringes the valid intellectual property rights of others. PendoTECH does not direct or instruct customers to use any PendoTECH product or system in combination with particular commercially available pre-sterilized manifolds, remotely operated pinch valves, as such use could possibly be covered by the claims of one or more patents. PendoTECH shall not be liable to BUYER in any manner with respect to Product sold. SELLER, PendoTECH MAKES NO REPRESENTATIONS OR WARRANTIES OF ANY TYPE, EXPRESS OR IMPLIED, AND EXPRESSLY DISCLAIMS AND EXCLUDES ANY REPRESENTATION OR WARRANTY OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE OR USE, NON-INFRINGEMENT OR WARRANTY ARISING FROM USAGE OF TRADE, COURSE OF DEALING OR PERFORMANCE."

For warranty information see our website at <http://www.pendotech.com/warranty>