




Pressure Sensor Cable Dust Cover/ Zero Simulator for PressureMAT®

The PendoTECH PressureMAT® is a process monitor that gives an indication of a high pressure situation when a PendoTECH Single Use Pressure Sensor™ is disconnected. This is to alert users of a potential un-intentional situation (versus remaining at zero) during a disconnect.

The Pressure Sensor Cable Dust Cover/Zero Simulator can be inserted into the PressureMAT monitor cables and the monitor will read zero. This prevents continually high pressure indication and alarm, particularly if the monitor is integrated to a central control system that logs alarms. In addition, this protects the gold pins when no sensors are connected, particularly if chemical cleaning solutions may be used around the PressureMAT cables.

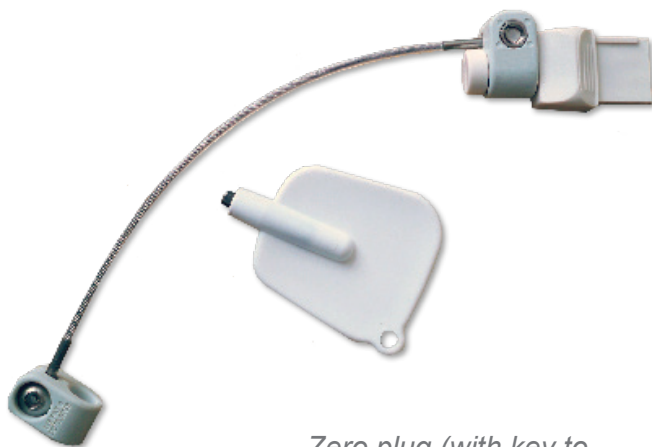

*Zero plug (with key to
attach to cable)*



*It comes with a tether so it remains
attached to the cable when not in use*



Ordering Information

PDKT-650-298CVR

Pressure sensor cable dust cover / zero simulator for PressureMAT

PendoTECH® and PressureMAT® are registered trademarks of PendoTECH, all rights reserved.

For warranty information see our website at <http://www.pendotech.com/warranty>